

1.0 MiCare Central Software Installer Manual

This document describes how to install the MiCare Central Software version 1.1.3.

1.1 Requirements

- Windows Server 2012 R2
- 1.4 GHz 64-bit processor minimum
- 512 MB RAM minimum
- 2 GB disk space minimum
- Internet access



Note: This installer:

- installs a new version of the MiCare Central Software,
OR
- upgrades the MiCare Central Software from version 1.1.1 to the latest version.

1.2 Before Installing

If you are upgrading from MiCare Central Software version 1.1.1 or higher, follow these instructions before installing.

1. Open Internet Information Services (IIS) Manager.
2. In Internet Information Services (IIS) Manager, right-click the MiCare website, then select **Remove**.


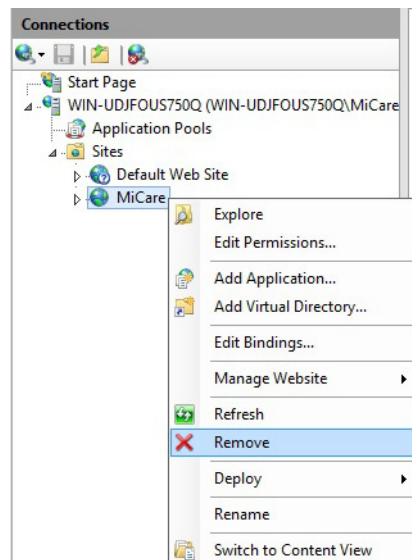


Figure 1 Internet Information Services (IIS) Manager

3. Go to section 1.3.

1.3 Installing


1. Double-click the file **MiCare Installer v1.1.3** to start the installation.
2. Follow instructions on the screen.
3. On the Email Server Configuration window, enter your own email server information, or leave the information as is. Consult your network administrator for assistance.

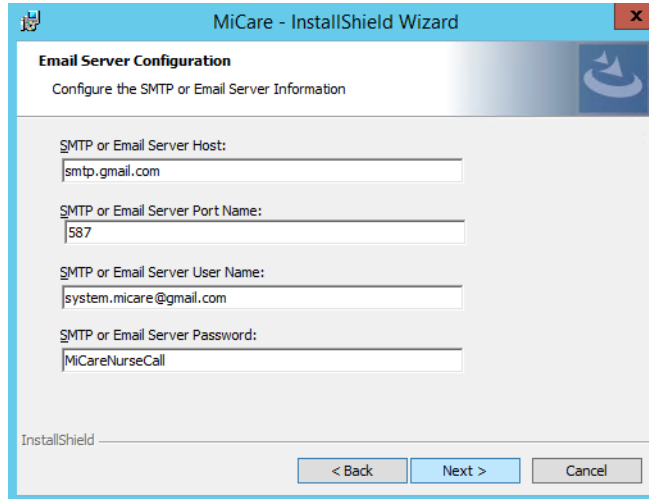

The screenshot shows the 'Email Server Configuration' window of the MiCare - InstallShield Wizard. The window title is 'MiCare - InstallShield Wizard'. The subtitle is 'Email Server Configuration' with the instruction 'Configure the SMTP or Email Server Information'. The window contains four text input fields: 'SMTP or Email Server Host' with the value 'smtp.gmail.com', 'SMTP or Email Server Port Name' with the value '587', 'SMTP or Email Server User Name' with the value 'system.micare@gmail.com', and 'SMTP or Email Server Password' with the value 'MiCareNurseCall'. At the bottom, there are three buttons: '< Back', 'Next >', and 'Cancel'. The 'Next >' button is highlighted in blue.

Figure 2 Email Server Configuration

4. On the Database Server window, click **Next**.

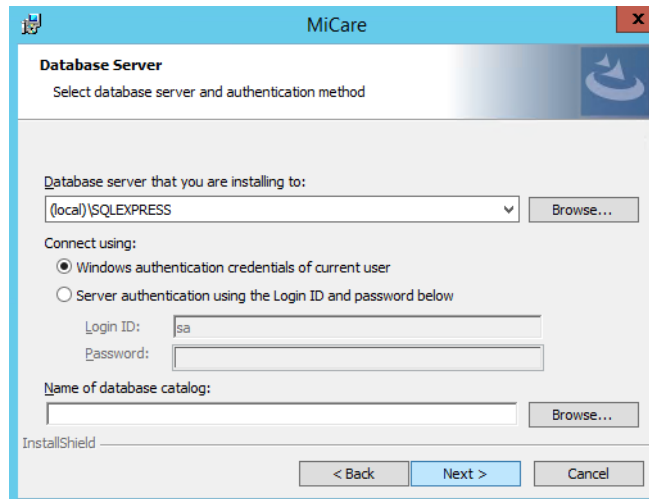

The screenshot shows the 'Database Server' window of the MiCare installer. The window title is 'MiCare'. The subtitle is 'Database Server' with the instruction 'Select database server and authentication method'. The window contains a dropdown menu for 'Database server that you are installing to:' with the value '(local)\SQLEXPRESS' and a 'Browse...' button. Below this, there are two radio buttons for 'Connect using:'. The first is 'Windows authentication credentials of current user' (selected) and the second is 'Server authentication using the Login ID and password below'. Under the second option, there are two text input fields: 'Login ID:' with the value 'sa' and 'Password:'. Below these, there is a text input field for 'Name of database catalog:' and a 'Browse...' button. At the bottom, there are three buttons: '< Back', 'Next >', and 'Cancel'. The 'Next >' button is highlighted in blue.

Figure 3 Database Server

5. When the menu Build MiCare Database GUI Tool appears (shown in Figure 4), follow the instructions below.
 - If you are installing the MiCare software on a new computer:
 - a. Press 1 then Enter.
 - b. Press 2 then Enter.
 - c. Press 3 then Enter.
 - d. Press 9 then Enter.
 - If you are upgrading the MiCare software:
 - a. Press 2 then Enter.
 - b. Press 3 then Enter.
 - c. Press 9 then Enter.

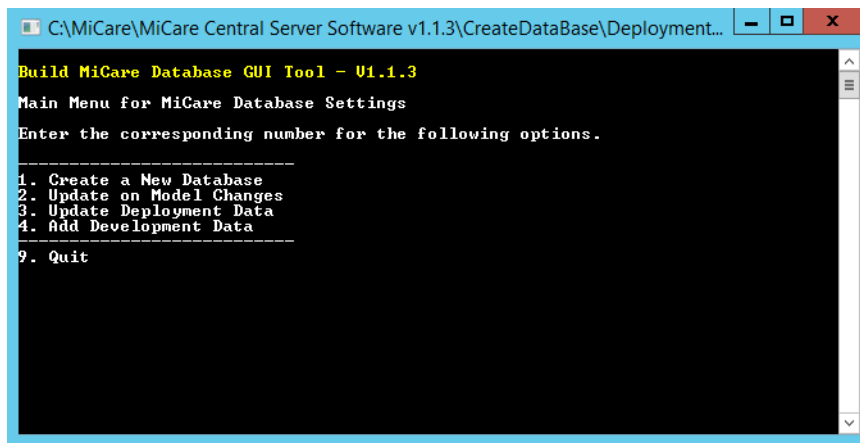


Figure 4 Build MiCare Database GUI Tool

6. Follow the instructions on the screen to complete the installation.

1.4 Contact Us

1.4.1 Canada and USA

Toll Free: 1-888-660-4655

Local: 905-660-4655

Fax: 905-660-4113

1.4.2 Website

<http://www.mircom.com>